

Alpha[®] Basic Liners

Product Highlights

- Designed specifically for low-activity use
- Flexible fabric comfortable for extended periods of knee flexion
- Pistoning minimized distally without a constrictive matrix
- Retrofits with common liner profiles


The Alpha Basic Liner was developed specifically for amputees with low activity levels. The liner has a Symmetrical 6 mm Gel Profile. WillowWood's Alpha Classic gel is paired with a flexible fabric to ensure easy donning and that the liners do not grab or pull on the knee when flexed. To minimize pistoning, the Alpha Basic Locking Liner utilizes the one-way stretch Select Fabric at the distal end of the liner without constricting a patient's residual limb.



- ✔ Transtibial
- ✔ Transfemoral

L-Code Suggestions: L5679 (cushion), L5973 (locking)

Liner Sizing:

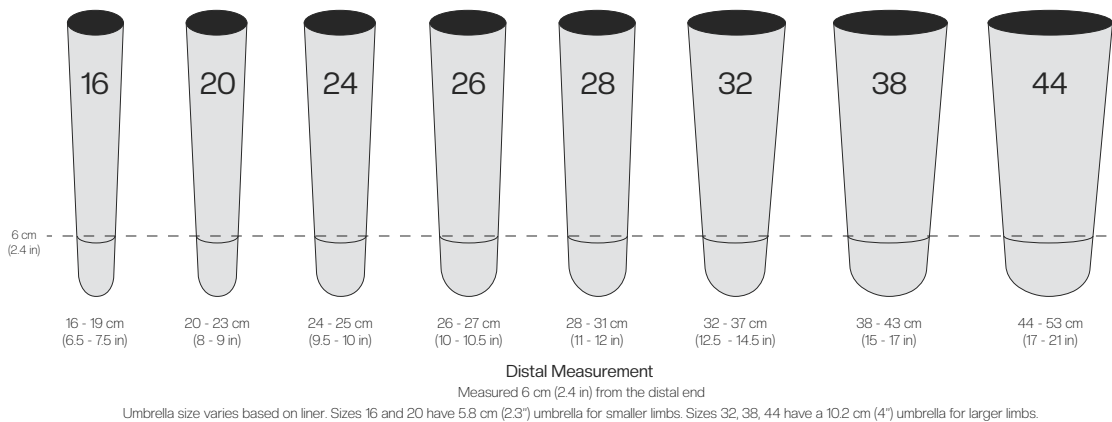

Part Numbers:

Size	Standard Locking Umbrella Item #	Large Locking Umbrella Item #	Cushion Item #
16	ABLL-16-6	-	ABCL-16-6
20	ABLL-20-6	-	ABCL-20-6
24	ABLL-24-6	-	ABCL-24-6
26	ABLL-26-6	-	ABCL-26-6
28	ABLL-28-6	-	ABCL-28-6
32	ABLL-32-6	ABKL-32-6	ABCL-32-6
38	ABLL-38-6	ABKL-38-6	ABCL-38-6
44	ABLL-44-6	ABKL-44-6	ABCL-44-6

Please Note: A drying stand is not included with the Alpha Basic Liner.