

Alpha[®] Hybrid Liners

Product Highlights

- Controls pistoning without a distal matrix
- Cushion liners ideal for use with elevated vacuum
- Size of accordion umbrella changes with liner size

The Alpha Hybrid Liner allows for easy bending and flexing of the knee, controlled pistoning, and easy donning. The liner's one-way stretch fabric effectively controls pistoning without a distal matrix and allows circumferential stretch while maintaining liner elongation within specified limits. A flexible panel incorporates an abrasion-resistant fabric allowing the limb to easily flex without pulling and tightening at the knee. This panel and the thin gel behind the knee allows for comfort while sitting or bending for extended periods of time. A new socket should be fabricated for best use of the Alpha Hybrid Liner with Progressive Gel style. The Alpha Hybrid Liner with Uniform Gel style is a solution for amputees currently wearing a 6 mm Uniform Alpha Classic Liner with Original Fabric and want to experience the benefits of this new liner.




✓ Transtibial


L-Code Suggestion: L5679 (cushion), L5973 (locking)

Liner Sizing:

Progressive Profile

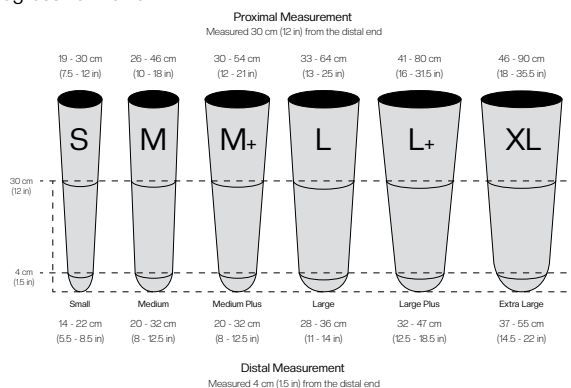

Uniform Profile

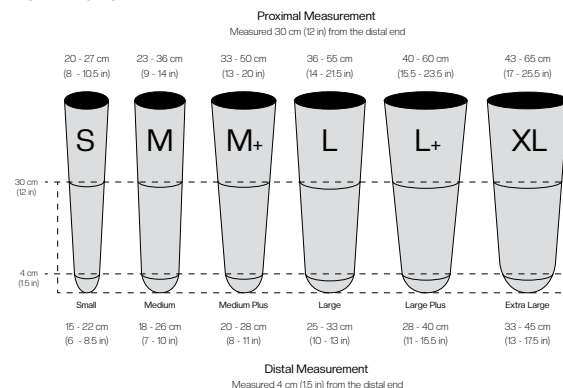

Part Numbers:

Size	Locking		Cushion	
	Progressive Item #	Uniform Item #	Progressive Item #	Uniform Item #
Small	H352-6393	H351-5363	H350-6393	H350-5363
Medium	H352-6390	H351-5360	H350-6390	H350-5360
Medium Plus	H352-6394	H351-5364	H350-6394	H350-5364
Large	H352-6396	H351-5366	H350-6396	H350-5366
Large Plus	H352-6398	H351-5368	H350-6398	H350-5368
Extra Large	H352-6397	H351-5367	H350-6397	H350-5367