

Alpha[®] Silicone

AK Liners

Product Highlights

- Non-greasy & non-tacky surface
- Controls pistoning without a distal matrix
- Accordion umbrella blends with a patient's anatomy


The Alpha Silicone AK Liner with Symmetrical Profile provides extra protection to the distal femur and allows the liner to be rotated for extended wear. The liner uses a proprietary, platinum-cured, medical grade silicone designed for both comfort and performance for transfemoral amputees. Our unique blend of silicone, Vitamin E, and skin conditioners results in a smooth surface that is non-greasy, non-tacky, and comfortable against the skin. The Alpha Silicone AK Liner uses Select Fabric that controls pistoning without using a distal matrix yet allows for circumferential stretch for easy donning.



✔ Transfemoral

L-Code Suggestion: L5679 (cushion), L5973 (locking)

Liner Sizing:

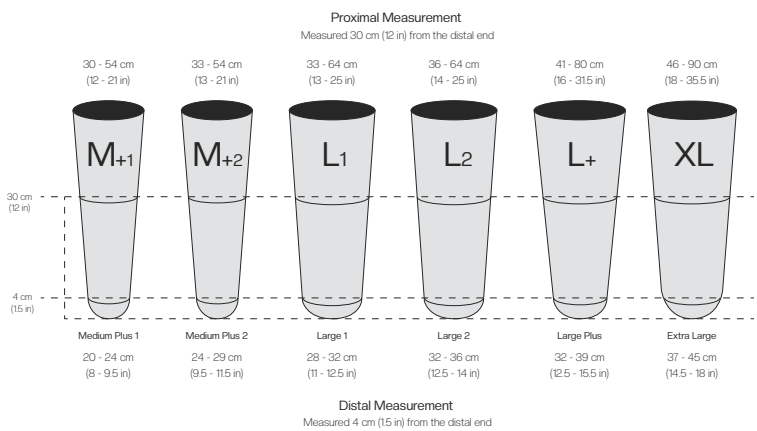

Part Numbers:

Size	Locking Item #	Cushion Item #
Medium Plus 1	S442-1694-1	S440-1694-1
Medium Plus 2	S442-1694-2	S440-1694-2
Large 1	S442-1696-1	S440-1696-1
Large 2	S442-1696-2	S440-1696-2
Large Plus	S442-1698-1	S440-1698-1
Extra Large	S442-1697-1	S440-1697-1