

Alpha[®] Silicone Liners

Product Highlights

- Non-greasy & non-tacky interface
- Controls pistoning without a distal matrix
- Size of accordion umbrella changes with liner size




WillowWood's Alpha Silicone Liner uses a proprietary, platinum-cured, medical grade silicone designed for comfort and performance. Our unique blend of silicone, Vitamin E, and skin conditioners results in a pleasing surface that is non-greasy, non-tacky, and comfortable against the skin. The Alpha Silicone Liner uses Select Fabric that controls pistoning without using a distal matrix. The Select Fabric allows for circumferential stretch making it easy to don for numerous limb shapes. The flexible knee panel provides enhanced knee flexion and a broad range of knee motion. The Alpha Silicone Liner retrofits with Alpha Hybrid Liners with a Progressive Profile.

✔ Transtibial

L-Code Suggestions: L5679 (cushion), L5973 (locking)

Liner Sizing:

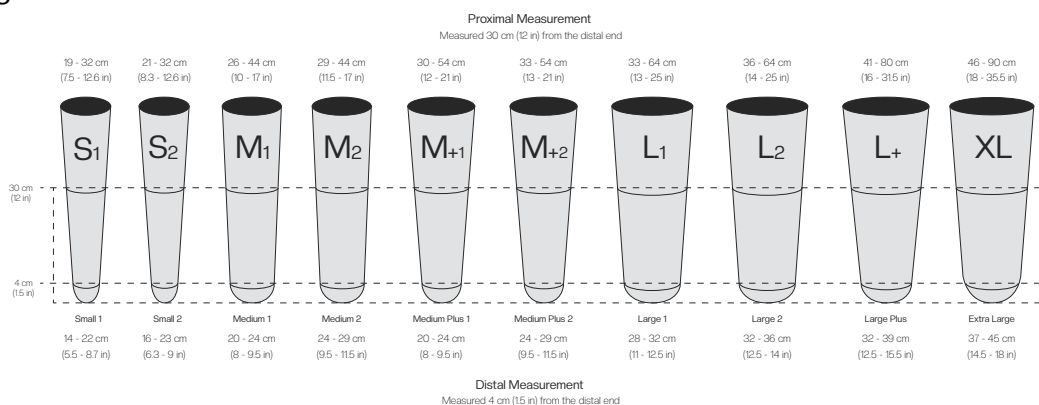

Part Numbers:

Size	Locking Item #	Cushion Item #
Small 1	S352-6593-1	S350-6593-1
Small 2	S352-6593-2	S350-6593-2
Medium 1	S352-6590-1	S350-6590-1
Medium 2	S352-6590-2	S350-6590-2
Medium Plus 1	S352-6594-1	S350-6594-1
Medium Plus 2	S352-6594-2	S350-6594-2
Large 1	S352-6596-1	S350-6596-1
Large 2	S352-6596-2	S350-6596-2
Large Plus	S352-6598-1	S350-6598-1
Extra Large	S352-6597-1	S350-6597-1