

Alpha[®] SmartTemp

Gel Liners

Product Highlights

- Proprietary blend of TPE gel and heat absorption technology
- Delays the onset of sweat
- Stabilizes skin temperature
- One-way stretch fabric controls pistoning

Sweat and heat within a prosthetic liner are long-standing complaints among amputees. Clinicians may now address these issues and provide their patients with the Alpha SmartTemp Gel Liner with WillowWood's proprietary blend of TPE gel and heat absorption technology. This technology allows the liner to absorb and store heat from a residual limb and delay the onset of sweat. The liner releases the stored heat as the body cools down, stabilizing the skin temperature to keep an amputee comfortable all day. With the reduction of sweat within the liner, chafing and skin breakdown may be reduced while suspension strength of the prosthesis may be improved.

The liner utilizes one-way stretch Select Fabric to control pistoning without the need of a distal matrix while allowing circumferential stretch for easy donning. The flexible knee panel of abrasion-resistant fabric permits comfortable, easy flexion without pulling and tightening at the knee. The panel and the thin gel behind the knee allow for comfort when sitting or bending for extended periods of time. Available in 6 mm Uniform and Progressive Profiles, the Alpha SmartTemp Gel Liner retrofits with other transtibial Alpha Liners of the same size, profile and distal end.




✔ Transtibial

L-Code Suggestions: L5679 (cushion), L5973 (locking)


Liner Sizing:

Progressive Profile

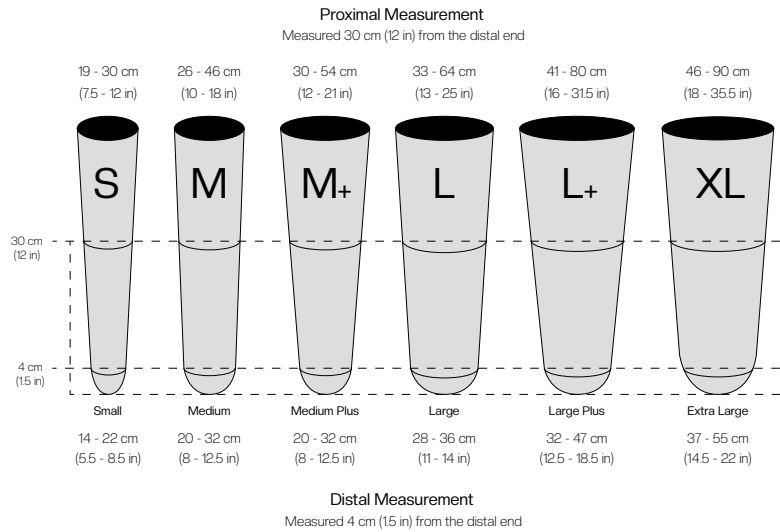

Patients wearing this liner may retrofit with other Alpha Liners of the same profile, base size, and umbrella style.

Uniform Profile

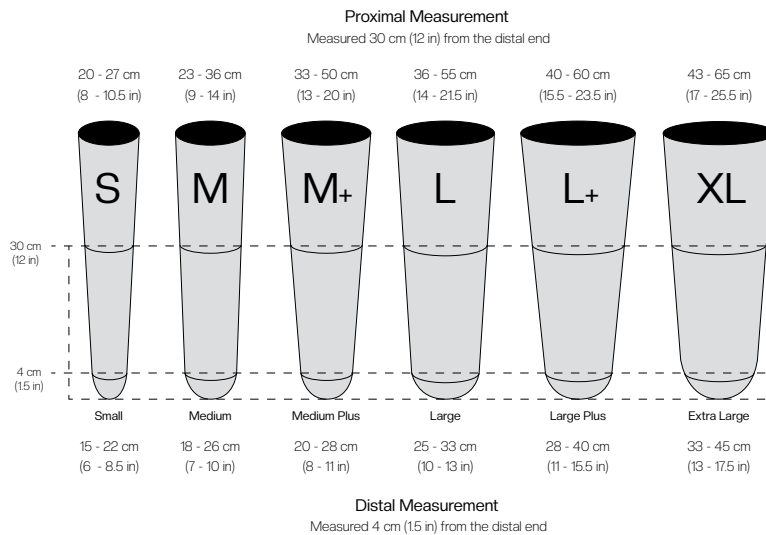

Patients wearing this liner may retrofit with other Alpha Liners of the same profile, base size, and umbrella style.

Part Numbers:

Size	Locking		Cushion	
	Progressive Item #	Uniform Item #	Progressive Item #	Uniform Item #
Small	P352-6393	P351-5363	P350-6393	P350-5363
Medium	P352-6390	P351-5360	P350-6390	P350-5360
Medium Plus	P352-6394	P351-5364	P350-6394	P350-5364
Large	P352-6396	P351-5366	P350-6396	P350-5366
Large Plus	P352-6398	P351-5368	P350-6398	P350-5368
Extra Large	P352-6397	P351-5367	P350-6397	P350-5367