

# Alpha® SmartTemp

## Silicone Liners

### Product Highlights

- Uses Outlast® technology - originally developed for NASA
- Delays the onset of sweat and reduces chafing
- Uses Select Fabric for a broad range of motion
- Available in Progressive and Uniform Profiles

The Alpha SmartTemp Silicone Liner featuring Outlast® combines silicone with the heat management technology of Outlast®. This revolutionary liner regulates heat by responding to fluctuations in skin temperature. The Alpha SmartTemp Silicone Liner absorbs and stores heat as it builds up in a residual limb delaying the onset of sweat before it starts. By releasing stored heat as the body cools, the liner stabilizes skin temperature to keep an amputee comfortable all day.

The Alpha SmartTemp Silicone Liner uses one-way stretch Select Fabric that allows circumferential stretch for easy donning while controlling pistoning. The liner is offered in both Progressive and Uniform Profiles which allows for retrofitting with Alpha Liners of the same profile, base size, and umbrella style.

The liner utilizes one-way stretch Select Fabric to control pistoning without the need of a distal matrix while allowing circumferential stretch for easy donning. The flexible knee panel of abrasion-resistant fabric permits comfortable, easy flexion without pulling and tightening at the knee. The panel and the thin gel behind the knee allow for comfort when sitting or bending for extended periods of time. Available in 6 mm Uniform and Progressive Profiles, the Alpha SmartTemp Gel Liner retrofits with other transtibial Alpha Liners of the same size, profile and distal end.




✓ Transtibial


L-Code Suggestions: L5679 (cushion), L5973 (locking)

## Liner Sizing:

### Progressive Profile

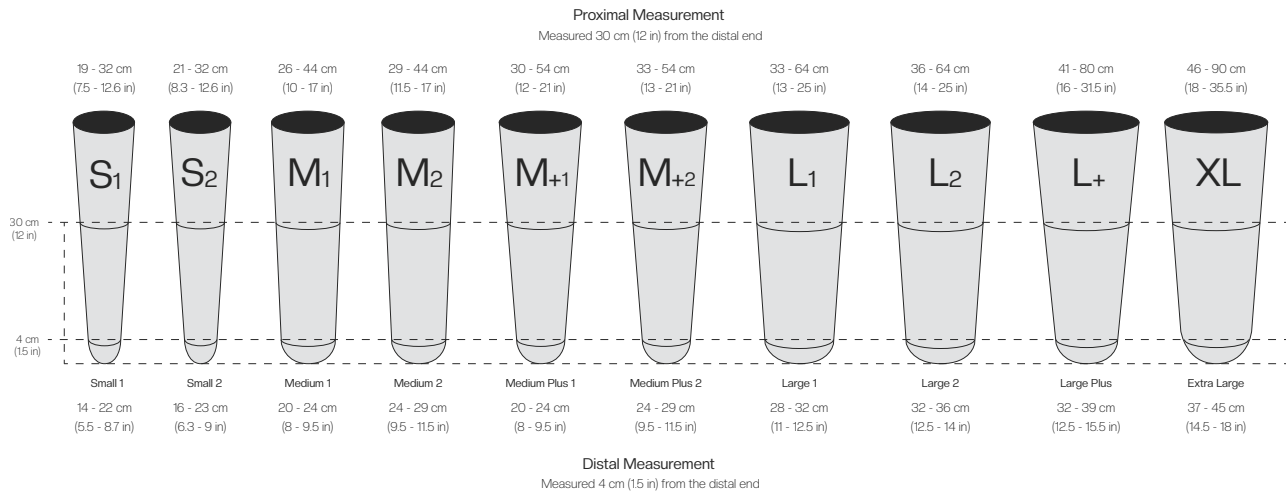

### Uniform Profile

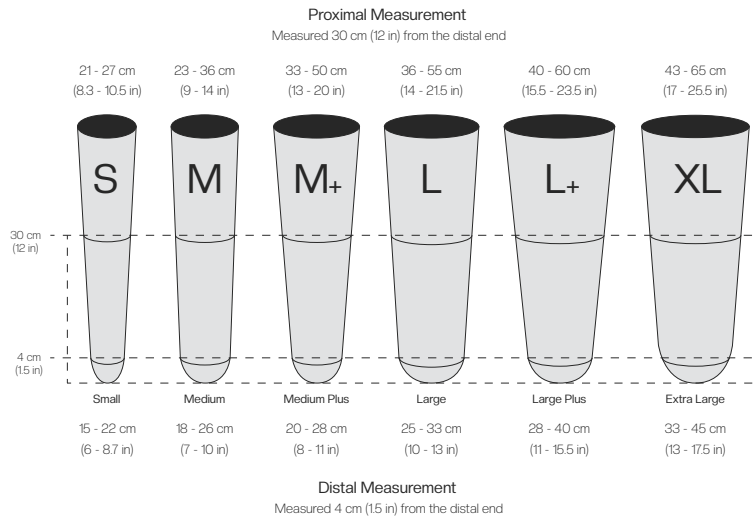

### Part Numbers:

| Size          | Locking Item # | Cushion Item # |
|---------------|----------------|----------------|
| Small 1       | T352-6593-1    | T350-6593-1    |
| Small 2       | T352-6593-2    | T350-6593-2    |
| Medium 1      | T352-6590-1    | T350-6590-1    |
| Medium 2      | T352-6590-2    | T350-6590-2    |
| Medium Plus 1 | T352-6594-1    | T350-6594-1    |
| Medium Plus 2 | T352-6594-2    | T350-6594-2    |
| Large 1       | T352-6596-1    | T350-6596-1    |
| Large 2       | T352-6596-2    | T350-6596-2    |
| Large Plus    | T352-6598-1    | T350-6598-1    |
| Extra Large   | T352-6597-1    | T350-6597-1    |

| Size/Profile  | Locking Item # | Cushion Item # |
|---------------|----------------|----------------|
| Uniform, 6 mm |                |                |
| Small         | T351-5363      | T350-5363      |
| Medium        | T351-5360      | T350-5360      |
| Medium Plus   | T351-5364      | T350-5364      |
| Large         | T351-5366      | T350-5366      |
| Large Plus    | T351-5368      | T350-5368      |
| Extra Large   | T351-5367      | T350-5367      |