

Alpha[®] TruSeal Liner

Product Highlights

- Medical grade, platinum-cured silicone
- Four integrated sealing fins
- Satin outer-finish, no lubricant needed
- Embedded fabric mesh
- Features a 12-month warranty

WillowWood's Alpha TruSeal Liner uses a proprietary, platinum-cured, medical grade silicone designed for comfort and performance. The liner has four integrated sealing fins designed to ensure confident suction suspension. With TruSeal's stepped fins and Volume Management Pads, a large range of fluctuation can be accommodated. The satin outer finish allows for easy donning, does not require messy lubricants and is easier to clean than a fabric covered liner. An inner topcoat is applied to reduce tack while controlling soft tissue compression and rotation. An embedded mesh fabric within the silicone resists abrasion and protects the liner from rips and tears. The Alpha TruSeal Liner is available in a Symmetrical Profile and two tapers (conical and cylindrical) to accommodate various limb shapes.




Alpha Silicone VMP, set of 2
Pad measures 9" x 4.5", 4.5 mm
& 7.5 mm thicknesses
ALA-VM-S2

✓ Transfemoral

L-Code Suggestion: L5679 (cushion)


Liner Sizing: Conical

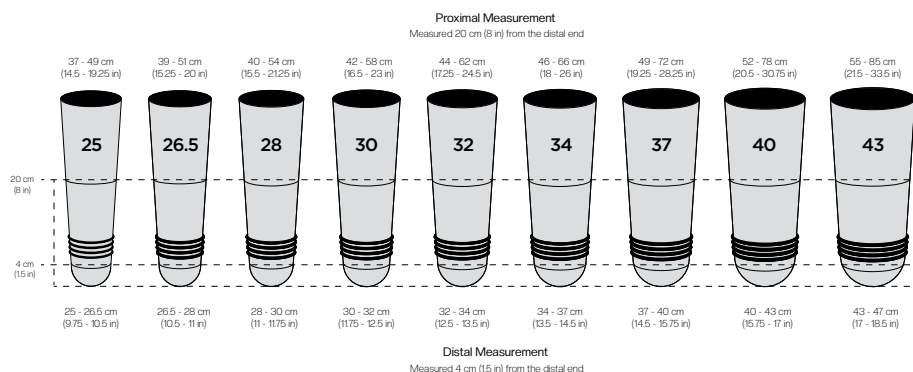

Part Numbers:

Size	Item #	Optional
25	S494-1125	ALA-VM-S2 Alpha Silicone Volume Management Pads, set of 2
26.5	S494-1126	
28	S494-1128	
30	S494-1130	
32	S494-1132	
34	S494-1134	
37	S494-1137	
40	S494-1140	
43	S494-1143	

Liner Sizing: Cylindrical

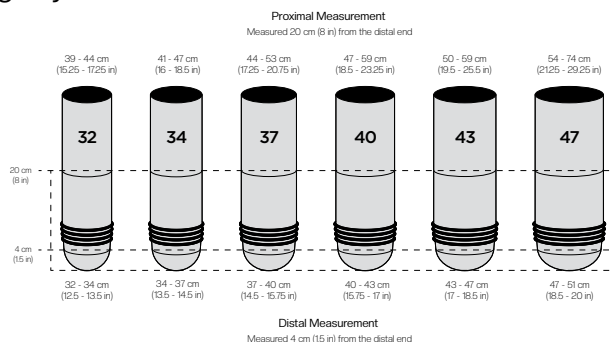

Part Numbers:

Size	Item #	Optional
32	S495-1132	ALA-VM-S2 Alpha Silicone Volume Management Pads, set of 2
34	S495-1134	
37	S495-1137	
40	S495-1140	
43	S495-1143	
47	S495-1147	