

# One Sleeve

## Product Highlights

- Pre-flexed at 22.5°
- Allows for easy knee flexion
- Use with WillowWood One® System for airtight seal

Designed for use with the WillowWood One® TT System, the One Sleeve is a pre-flexed internal sleeve that allows for a comfortable fit without limiting the knee's range of motion. The fabric is adhered to the proximal and distal areas of the sleeve but left loose in between to allow for knee flexion without creating friction on the knee. The fabric-free proximal portion of the sleeve may be trimmed to suit patient needs. The sleeve is designed not to bunch behind the knee when standing, sitting or bending.

Unique among WillowWood's sleeve products, the 'fins' on the One Sleeve establish an airtight seal between the flexible inner socket and the laminated socket when used as part of the WillowWood One® TT System.


Warranty: 3 months from date of patient fitting



✔ Transtibial

L-Code Suggestion: L5685

### Sleeve Sizing:

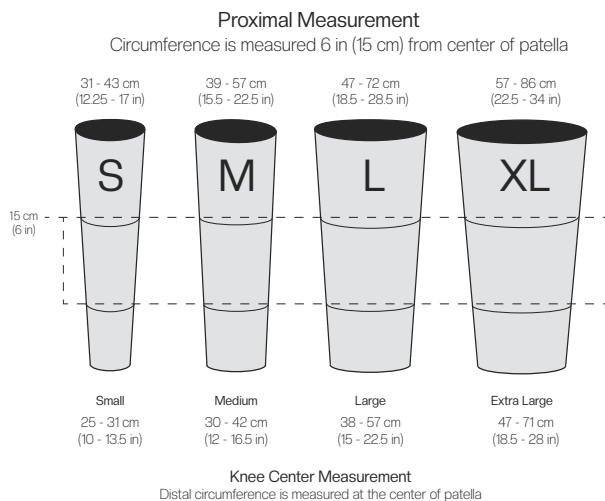

### Part Numbers:

Product	Item #
One Sleeve, Small	ONE-SLEEVE-TT-SM
One Sleeve, Medium	ONE-SLEEVE-TT-MD
One Sleeve, Large	ONE-SLEEVE-TT-LG
One Sleeve, Extra Large	ONE-SLEEVE-TT-XL